


Serum Creatinine: IDMS vs Legacy Calibration

Central Maine Medical Center

Vitros Systems

Sep-2007 with precision to 2 decimals

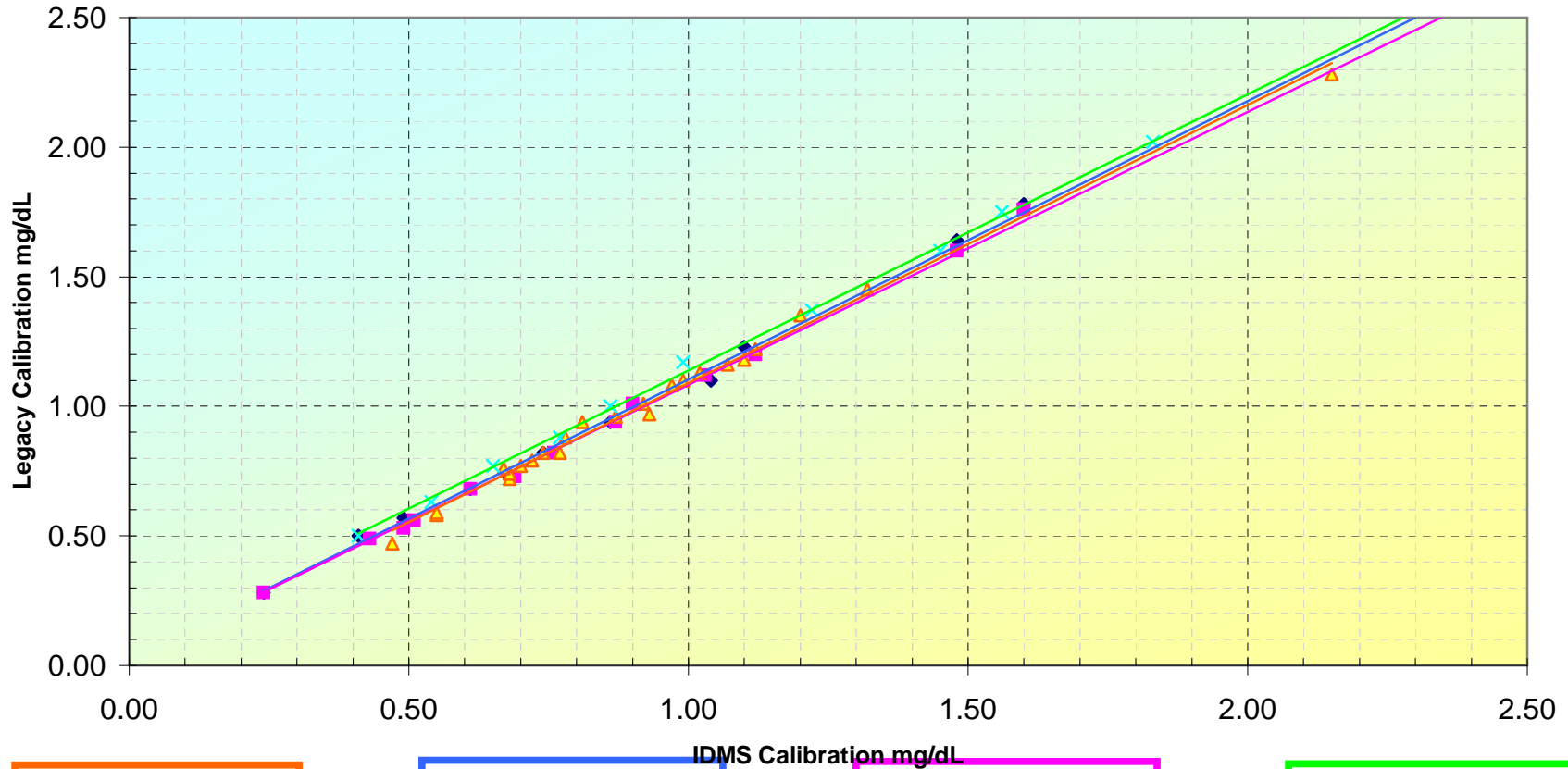

Rumford Vitros 250
 $y = 1.0728x + 0.017$
 $R^2 = 0.9945$

CMMC Fusion 007
 $y = 1.074x + 0.03$
 $R^2 = 0.9992$

CMMC Fusion 005
 $y = 1.0523x + 0.0313$
 $R^2 = 0.9994$

BH Vitros 350
 $y = 1.0653x + 0.0732$
 $R^2 = 0.9994$